









OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.





#Letsgivemore♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.





Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or.air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® Centrifuges

Blood Washing Centrifuge

Model YR142-2



Product Description:

- ✓ Microcomputer control, DC brushless motor, stable operation, low noise and high speed accuracy.
- ✓ Touch panel, programmable operation, operation data can be set according to demands and stored automatically.
- ✓ LED display. user-friendly interface, simple and convenient operations.



Technical Specifications:

Model	YR142-2		
Rotor type	HLA	SERO	
Rotor capacity	(0.25ml-1ml)x12	7mlx12	
Max RCF	2000xg	1000xg	
Drive	Micocomputer control, brushless DC motor drive		
Time range	Program fixed		
Noise	<60db(A)		
Power supply	AC110V/AC220V, 50Hz-60Hz		
Consumption	150w		
Dimension(L×W×H)	330mm×420mm×280mm		
Packing size(L×W×H)	430mm×520mm×390mm		
Net weight (No rotor)	15kg		









KALSTEIN® Centrifuges

Blood Washing Centrifuge

Model YR142-2




Rotor	Key	RCF Fixed	Time Fixed	Effective Solution
HLA ROTOR	1	2000xg	180s	Lymphocyte separation incubated cell separation
	2	1000xg	3s	Hematob last exenterate(thrombin disposal)
	3	1000xg	60s	Lymphocyte Washing
SERO ROTOR	1	500×g	60s	Blood grouping, image of hemagglutination reaction
	2	1000xg	15s	Cross matching, and the testing of antistaphyoloysin
	3	1000xg	60s	Hematocrit washing, the extraction of serum and plasma







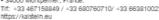









#LETSGIVEMORE

WITH THE ADCQUISITION OF A KALSTEIN EQUIPMENT

YOU MAKE A CONTRIBUTION TO:





https://onetreeplanted.org



biomédicale

https://www.humatem.org/



https://www.maniapurefoundation.org/




ANY DOUBT?

Contact us!

Articles Kalstein

(Social Networks and Official Websites)



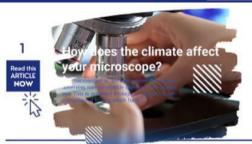



Article of the week: How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
https://kalstein.eu/articles





Article of the week: What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites https://kalstein.eu/articles



PARIS - FRANCE HEADQUARTER

2 Rue Jean Lantier 75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89 **Cel**: +33 (0) 6 80 76 07 10

sales@kalstein.eu

https://kalstein.eu/ https://www.kalstein.fr/


Kalstein Training Videos

(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...

23 views • 3 years ago

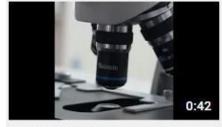

Kalstein's Microscopes

445 views • 1 week ago

