




LABORATORY LINES

DRYING OVEN



Kalstein®

A different accompaniment, at your service

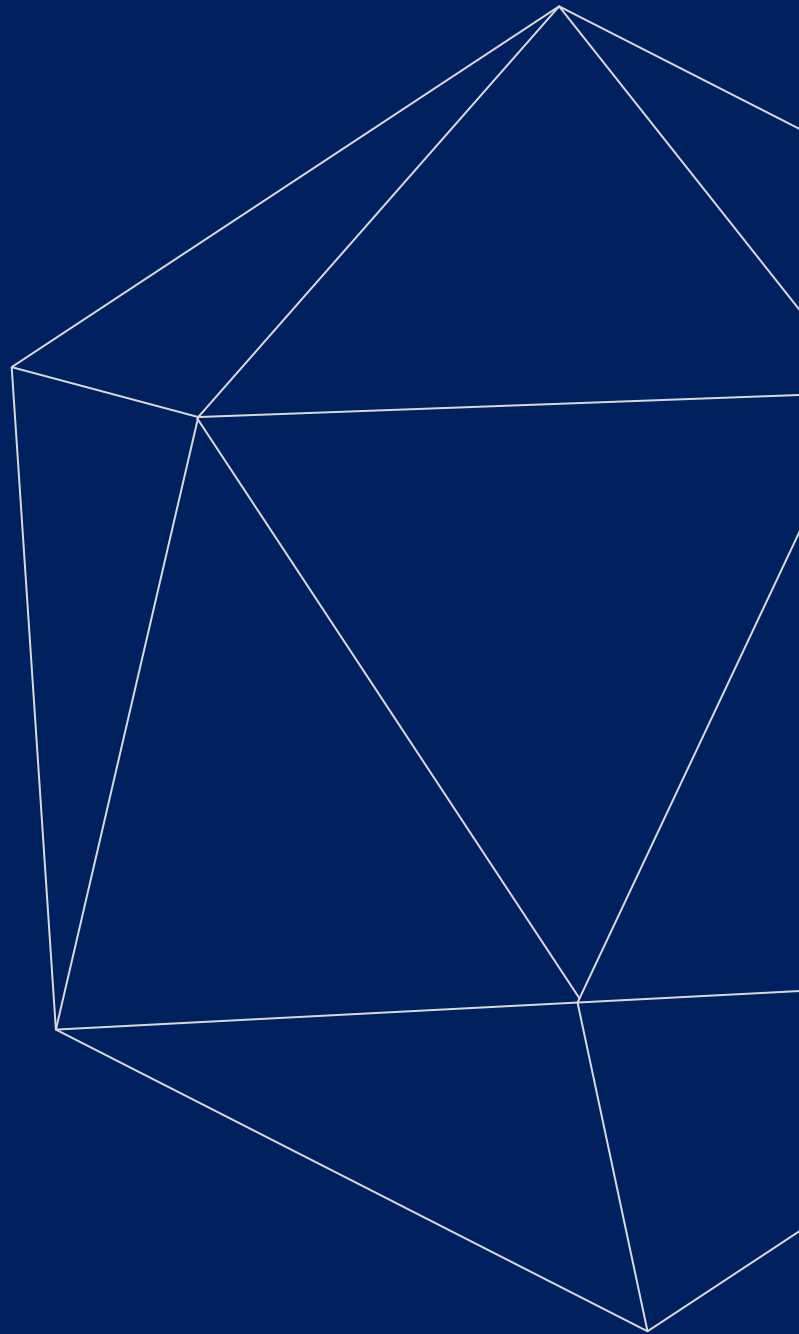



Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore ♥

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



Drying Oven

Serie YR101-5AS (1000L)



Two kinds of temperature controller; pointer controller and Microcomputer intelligent controller type; a forced convection, low noise hot air circulation system; shelf space and shelf number can be adjustable; double-layer tempered glass observation window.

Suitable for industrial factories, labs, research and other departments for drying, baking, heat treatment heat, it is not for volatile, flammable, explosive and other items. In order to avoid causing an explosion.

Suitable for routine testing of industrial and mining enterprises, basic laboratory and workshop is not recommended for a long-time continuous use. The requirement of higher temperature uniformity and large water content of the sample is not recommended.

Product introduction

Provide for the industrial factory, labs, research and other departments for Drying, baking, melting and sterilizer use.

Features :

- ✓ Three kinds of controller types are: Pointer controller, LED microcomputer controller and Microcomputer intelligent type.
- ✓ The shell is made of high-quality cold-rolled steel with metallic paint spray. Nice overall colour design suitable for all research labs and workshop production.
- ✓ With double-layer glass observation window, it is easy to watch the materials in the .
- ✓ Two models material inner chamber mirror stainless steel and ordinary cold-rolled stainless steel.
- ✓ Air circulation system consists of low noise fans and an air cycle duct to make the temperature more uniform. The reasonable air circulation ensures the temperature uniform in the inner chamber and fast warm-up time.
- ✓ High temperature resistant silicone for sealing the door, high-temperature heat pipes for the heating element to ensure safety and longer life.
- ✓ Fan forced convection for YR A series oven air circulation.



#Let'sgive more ♥

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



Drying Oven

Serie YR101-5AS (1000L)

 Energy-saving and environmental friendly



Model	Voltage (V)	Temp range (°C)	Accuracy (°C)	Fluctuation (%)	Power (KW)	Liter	Inner size (mm) (H·W·D)	Overall dimension (mm) (H·W·D)	Packing size (mm) (H·W·D)	Net weight (KG)	Shipping weight (KG)
YR101-0(AS)	220	RT+5 ~300	±1	±1	0.8~1.6	42	350*350*350	590*660*520	710*740*600	34.8	37.2
YR101-1(AS)	220			±1	0.8~1.6	70	450*450*350	690*760*520	810*840*600	41.5	45
YR101-2(AS)	220			±1	1.6~2.4	136	550*550*450	790*870*640	910*950*720	60	64.5
YR101-3(AS)	220			±1	2~3	225	750*600*500	1050*980*700	1270*1060*780	80	84.5
YR101-4(AS)	380			±2	3~4.8	640	1000*800*800	1300*1210*1000	1420*1290*1080	162	168
YR101-5(AS)	380			±2	4~6	1000	1000*1000*1000	1300*1420*1200	1420*1500*1280	189	197



#Letsgivemore ♥

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



DRYING OVEN



Temperature display controller



Vent hole



Cooling motor



Stainless steel inner



Double-layer tempered glass



#Letsgivemore ♥

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S




DRYING OVEN




Stainless steel inner (type with "s")



Zinc-painted steel inner

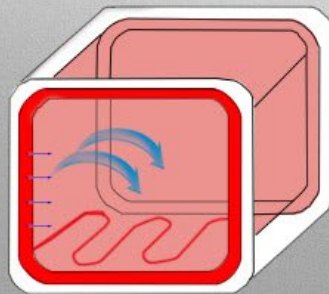

⊕ Airflow circulation system

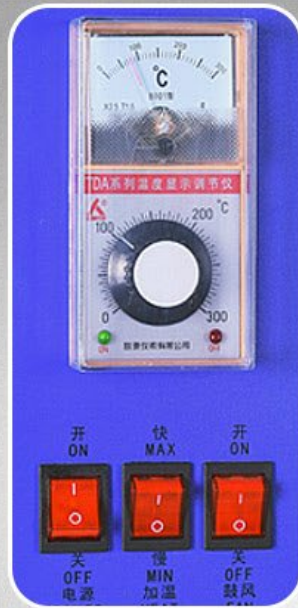

Pointer temperature control



Digital temperature control



#Letsgivemore

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.




All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



Optional accessories



1. Intelligent temperature controller

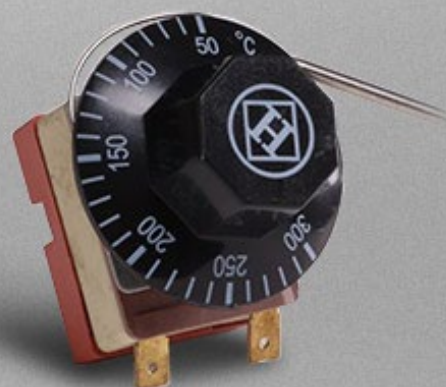

2. Independent temperature controller

Back



#Letsgivingore

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tit: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



DRYING OVEN

Scope



Jade processing



Lab glassware drying



Medical use



Auto headlight remould



Herbs drying



Circuit board aging test



#Letsgivemore

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



Articles Kalstein
(Social Networks and Official Websites)



Kalstein Blog
Not just be informed, be fully informed
The 3 articles of the week!
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



1
Read this ARTICLE NOW
How does the climate affect your microscope?
A microscope is an instrument that allows observing and measuring objects that are too small to be seen with the naked eye. This is achieved through a complex system of lenses and mirrors, which form...

Kalstein Blog
Article of the week:
How does the climate affect your microscope?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW
What is the importance of the laboratory refrigerator?
A refrigerator is a device that is one of the most common pieces of laboratory equipment. It is used to maintain, in a controlled environment (refrigerated space), various fluids and substances...

Kalstein Blog
Article of the week:
What is the importance of the laboratory refrigerator?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...

23 views • 3 years ago



Kalstein's Microscopes

445 views • 1 week ago



MONTPELLIER - FRANCE
HEADQUARTER
Optimum Business Center,
450 Rue Baden Powell,
34000 Montpellier - France.

Fax: +33 (0) 4 67 15 88 49
Tif: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

Official Websites for Latin America.

